

FOB Gold®  
LINE

# SENTiFOB®

Automated system for  
Fecal Immunochemical Testing

**SENTINEL**  
DIAGNOSTICS

# SENT*i*FOB®

## Easy, safe and automated Fecal Immunochem



### PIERCING FUNCTION

for FIT testing

- Automatic **PIERCING** function
- Sample probe with **LEVEL SENSOR**
- Patented and dedicated **PIERCEABLE TUBES**
- **READY AND EASY** to use for fecal sampling
- Tubes openable only at the **COLLECTION** side
- **SAFE AND HYGIENIC** for the operator
- **READY FOR THE ASSAY:**  
pierceTubes are loaded directly  
on the instrument for automatic analysis



### SAMPLES AND REAGENTS management

- On-board total samples: **72 TUBES**
- Sample loading capacity:  
**3 SAMPLE RACKS** of 24 each
- **CONTINUOUS LOADING** of samples
- Reagent loading capacity:  
**MULTIPLE VIAL POSITIONS**
- On board total tests:  
**480 DETERMINATIONS**
- External **BARCODE READER**  
for sample identification



# ical Testing



## **FAST AND FULLY AUTOMATED** analyzer

- **WALK-AWAY** system, from start to results
- Throughput up to **100 TESTS/H**
- First result after only **14 MIN**
- **IMMUNOTURBIDIMETRIC** method
- **QUANTITATIVE** determination of FOB/FIT
- **AUTOMATIC DILUTION** of calibrator



## **EASY AND INTUITIVE** software interface

- **USER-FRIENDLY AND COMPLETE** software
- **EASY AND INTUITIVE** interface on real monitor
- **VISUALIZATION** of reaction curve
- **COLOR IDENTIFICATION** of sample results
- **COMPLETE** QC management
- **GUIDED AND SIMPLE** maintenance menu
- **EASY EXPORT** and **PRINTING** of results
- **LIS** bidirectional connectivity



**FOB Gold®  
LINE**

**SENTiFOB®**



Ref. 1156188

**SENTiFIT® pierceTube**

100 pcs

- Buffer formulation with enhanced hemoglobin stability



Ref. 1156009

**FOB Gold® NG Wide**

R1: 1 x 20 mL

R2: 1 x 20 mL

- Ready to use
- Open vials stability: 30 days



Ref. 1157501

**FOB Gold® Screen Diluent**

2 x 20 mL

- Ready to use
- Open vials stability: 30 days



Ref. 1157103

**FOB Gold® Control Set Wide**

2 x (1 x 2) mL

- Ready to use
- Negative and positive level
- Open vials stability: 30 days



Ref. 1157401

**FOB Gold® Calibrator Wide**

Cal 1 x 2 mL + Dil 1 x 5 mL


- Ready to use
- Calibration stability: 30 days

**Sentinel CH. SpA**

Via Robert Koch, 2  
20152 Milan - Italy  
Phone: +39 02 3455141  
Fax: +39 02 34551464  
sentinel@sentinel.it  
www.sentinel diagnostics.com



**SENTINEL**  
DIAGNOSTICS

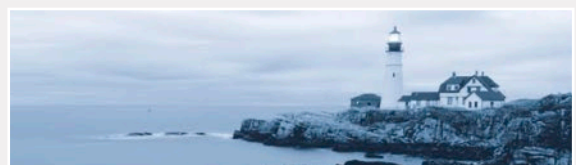

WATCHING OVER LIFE